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
Innovative design with four heating massage heads can produce kneading massage effect. It can deeply massage fatigue muscle tissue. Accelerate circulation of blood. Ease the pain, relax mood, and make you more comfortable to fully enjoy the fun of massage.

Product Parts and Accessories

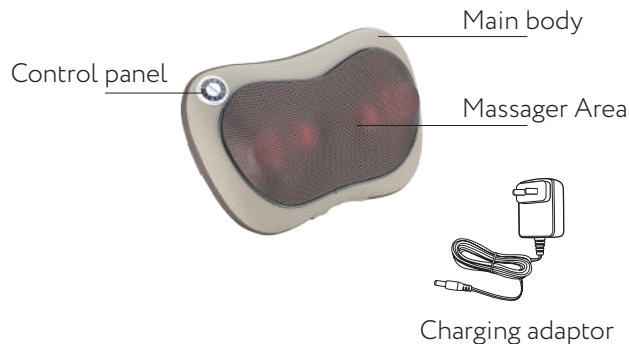Product Features

- Suitable for whole body massage, especially for the back, calf and the neck.
- With 4 massage heads with clockwise and anti-clockwise rotation kneading mode provides the comfortable massage experience.
- Automatically power off after 15 mins use.
- Wireless style with lithium battery: 2000mAh, power 12W.
- Adjustable massage speed : slow and fast
- Accessories: a charge adaptor

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Neck

Operating Instructions

- Rating massage time: 15mins, automatically power off after 15mins use; Press the On/Off button again if want to continue.

ON/OFF Button

- Press the On/Off button the first time, start massage with heating function.
- Press the On/Off button the second time, reverse rotation.
- Press the On/Off button the third time, turn off the heating function but keep massage.
- Press the On/Off button the fourth time, turn off the device, stop massage.

Strength adjustment Button

- Stay at the strong strength when turn on the massager, press the button can turn to weaker strength.

Charging Instruction

- Lithium battery inside the massager, full charge before first using.
- The indicator on the adaptor will be red when charging, and it will turn to be green when full charged. Unplug the adaptor to use the massager.
- In daily use, please charge it when the massager can't turn on or the massage strength turns weak.

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Caution

1. Unplug the device from the socket while not using.
2. Do not stand on the device.
3. The massager has a hot compressing function, and temperature can reach at about 45 C, users not sensitive to change in temperature should be careful in use.
4. Don't use in the bathroom or high humidity of the place; don't immerse in water or near the water, otherwise it will cause a short circuit or electric shock.
5. Do not use the device on the heating article, such as Electric blanket.
6. Stop using the device while the power cord or other element is damaged.
7. Do not disassemble the device.
8. Do not use the device in sleep and after drunk.
9. Children and the disabled using the device should be under the supervision.

Warning

The following conditions must consult a doctor before use; patients with heart disease, tumor patients, pregnant women and physiological period, spine malformation or abnormal person, have been treated in department of history with osteoporosis, implantation of cardiac pacemaker patients are prohibited the use of massage.

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Troubleshooting and technical performance index

If there is any problem when you use the device, please follow the instructions to check by yourself first. If still can't solve the problem, please cut off the power and contact the local service center for help. Do not disassemble or fix the device by yourselves.

Symptom	Possible Cause	Solution
The product is not working	Please check if you have press the power on switch button	Turn on the power switch
	Check if the battery has been run out	Charging using the original adaptor
	Check if time to the default setting time (15 mins)	Restart the power switch
	Weighing too much on the device	Unplug the device and remove the overload things
Device over heating	Long time use	Stop using the device for 2 hours
Power cord damaged, smoke comes out from the device		Stop using and contact with the local service center

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